

Swing™ Parts Guide

- | | |
|------------------------------|---|
| A. Motor unit | H. Shoulder/neck strap |
| B. Pump body | I. Power cord |
| C. PersonalFit™ breastshield | J. 2 - 5 oz (150 mL) BPA-free
breastmilk bottles |
| D. Tubing | K. Bottle stand |
| E. Valve | L. 2 bottle caps |
| F. Membrane | M. Travel cap |
| G. Drawstring storage bag | N. Slow flow wide base nipple |



medela 

www.medela.com

Medela, Inc., P.O. Box 660, 1101 Corporate Drive, McHenry, IL 60051-0660
Phone: (800) 435-8316 or (815) 363-1166 Fax: (815) 363-1246 Email: customer.service@medela.com

Swing, PersonalFit and #1 Choice of Hospitals & Mothers are trademarks and
2-Phase Expression and Medela are registered trademarks of Medela, Inc.

1908205 A 1210


© 2010 Medela, Inc.

medela 



Scanned + Archived
By:
FindMyManual.org

Swing™ Quickstart

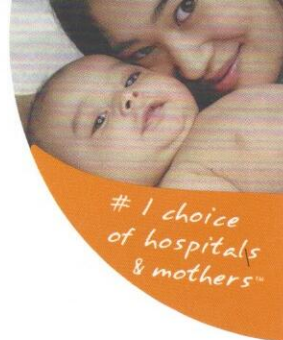


Swing™ QuickStart Instructions


See Medela's Swing instruction booklet for detailed information on assembly, cleaning and pumping.




For use with a Swing Breastpump

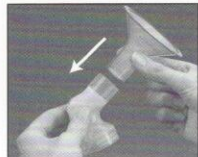


- 1 Make sure breastshield and connector are securely attached.

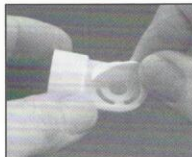


- 2 Snap membrane onto valve until membrane lies completely flat.

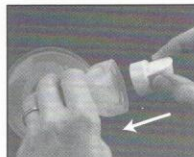

- 3 Push assembled valve & membrane onto connector.

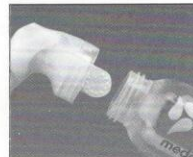

- 4 Screw bottle into connector.

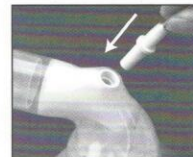


- 5 Insert round, light yellow end of tubing into opening on the pump body.



- 6 Insert the yellow, TRIANGLE end of tubing into the TRIANGLE port onto the motor.


NOTE: Be sure not to insert triangle end of tubing into the AC adapter port - Swing will not turn on if this is done. The AC adapter port is round in shape.

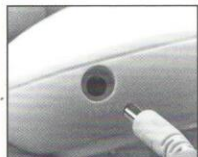


- 7 Plug the metal end of the power cord into the AC adapter port as shown. Plug the cord into a standard electrical outlet.




- 8 Center the breastshield over your nipple.
NOTE: Do not tilt or overfill containers while pumping. Take care not to kink tubing while pumping.

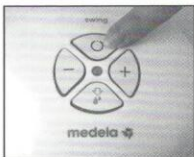

- 9 Press the on/off button to start Swing in the Stimulation Phase. **NOTE:** Do not hit the on/off button quickly or hold the on/off button down longer than one second or Swing will not turn on. Simply press the on/off button once slowly and release.

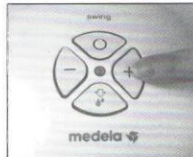

- 10 Adjust the vacuum to your comfort level while pumping by pressing the increase (+) or decrease (-) button. Holding down these buttons will adjust vacuum more rapidly. (The green LCD light will flash rapidly when you have reached minimum and maximum pumping settings.)

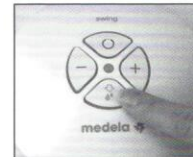

- 11 After two minutes, the Stimulation Phase will automatically change to the Expression Phase.



Properly connected tubing and AC adapter cord.